




# Liquid Biopsy

Accelerate Your Circulating Biomarker Research

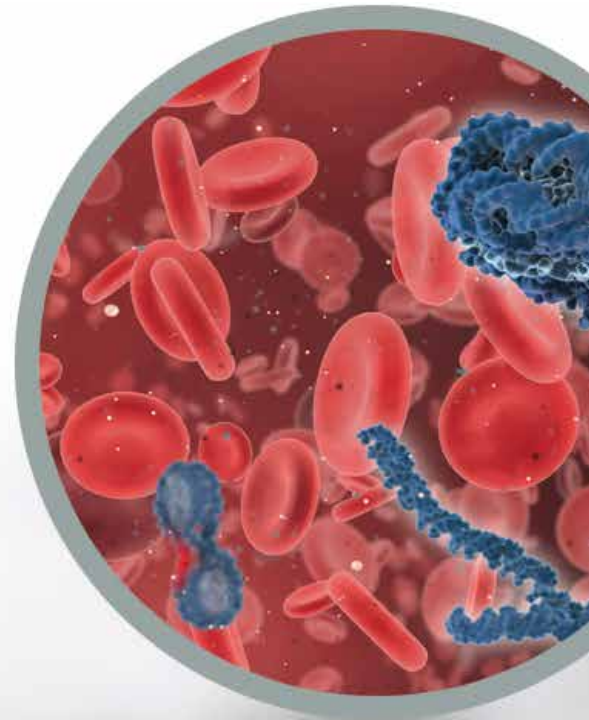

Sample to Insight

# Cross the bridge from discovery to insight

Liquid biopsy is a new non-invasive approach to help you detect disease biomolecules in blood, urine and other body fluids. It is a useful method when:

- Your tissue sample is limited
- Your sample contains insufficient tumor tissue
- Tumors are hard to reach
- Regular monitoring is needed


Liquid biopsy begins with sample processing and demands high efficiency for successful downstream detection. To unlock the molecular biomarkers in body fluids, QIAGEN® has sample processing solutions for all major liquid biopsy approaches – including automation.

A workflow with dedicated NGS library prep solutions ensures optimized conversion rates for either whole genome or targeted DNA sequencing and allows the detection of low-frequency variants with high confidence.

QIAGEN's Bioinformatics Advanced Testing Solution combines state-of-the-art data analysis and interpretation. Optimized, streamlined end-to-end workflows enable high sensitivity in detection, filtering and interpretation of variants.

ccfDNA/RNA = circulating cell-free DNA or RNA.

## Blood collection and stabilization

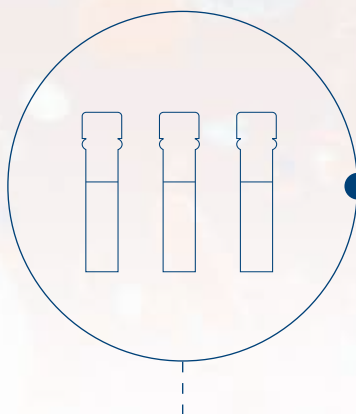


- Ensure long-term stability of markers
- Eliminate analysis errors caused by transportation and storage

### ccfDNA

- PAXgene® Blood ccfDNA Tubes

## Sample preparation

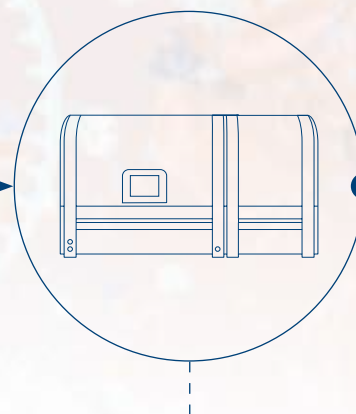

- Explore the biomarker universe with circulating nucleic acids, exosomes and CTCs
- Achieve highest yields and purity from precious samples

### ccfDNA

- QIAamp® MinElute® ccfDNA Kits
- EZ1® ccfDNA Kits
- QIASymphony® PAXgene Blood ccfDNA Kit
- QIAcube®
- EZ1 Advanced XL
- QIASymphony

### ccfRNA

- miRNeasy Serum/Plasma Advanced Kit
- exoRNeasy Serum/Plasma Kit

### ccfRNA & DNA (total nucleic acids)

- QIAamp ccfDNA/RNA Kit


### Exosomes

- exoEasy Maxi Kit
- exoRNeasy Serum/Plasma Kit
- miRCURY Exosome Kits

### Circulating tumor cells

- AdnaTest CancerSelect

## NGS analysis

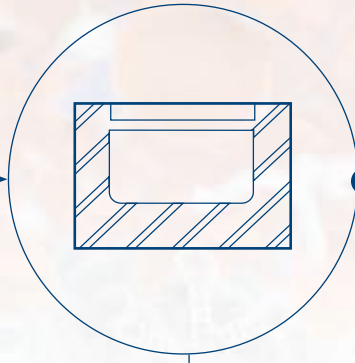

- Focus on relevant markers with targeted sequencing
- Discover the biomarker space with best-of-class platform-agnostic NGS solutions

### DNA sequencing

- QIAseq™ Targeted DNA Panels
- QIAseq cfDNA Library Kits

### miRNA sequencing

- QIAseq miRNA Library Kits

## Data analysis



- Identify variants and determine changes in expression with high accuracy
- Validate results with integrated visualization

- Biomedical Genomics Workbench
- CLC Genomics Workbench


## Interpretation



- Generate biological hypotheses
- Identify and prioritize variants for follow-up

- Ingenuity® Variant Analysis
- Ingenuity Pathway Analysis (IPA®)

## PCR-based analysis

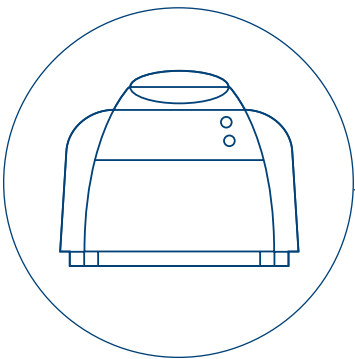

- Identify disease-relevant gene expression signatures
- Achieve meaningful results with wet bench-validated assays

### RNA biomarker profiling

- miRCURY PCR Panels
- RT<sup>2</sup> Profiler Arrays
- RT<sup>2</sup> lncRNA Profiler Arrays
- QuantiNova® PCR Kits

### Circulating tumor cells

- AdnaTest CancerDetect

## Collection and stabilization

# Ensure accurate preservation

## PAXgene Blood ccfDNA Tubes (RUO)


Avoid variation in the quality and integrity of circulating, cell-free DNA (ccfDNA) isolated from blood specimens. PAXgene Blood ccfDNA Tubes contain a reagent that prevents gDNA release to effectively stabilize ccfDNA levels in the plasma without impacting downstream assays. Pair seamlessly with the QIAamp MinElute ccfDNA and EZ1 ccfDNA kits, or benefit from optimized system performance and primary tube handling options when using the QIA Symphony PAXgene Blood ccfDNA Kit.

Novel noncrosslinking stabilization reagent compatible with a wide spectrum of downstream applications.


- Accurately detect and quantify ccfDNA from blood samples stored in tubes for up to 7 days at room temperature (15–25°C)
- Test for methylated circulating tumor DNA even in challenging assays

Unparalleled stabilization of red blood cells and minimized hemolysis during whole blood storage.

- Maximize plasma recovery
- Mitigate the risk of gDNA background
- Eliminate a second centrifugation step in plasma prep
- Allow possibilities for gDNA analysis from the buffy coat

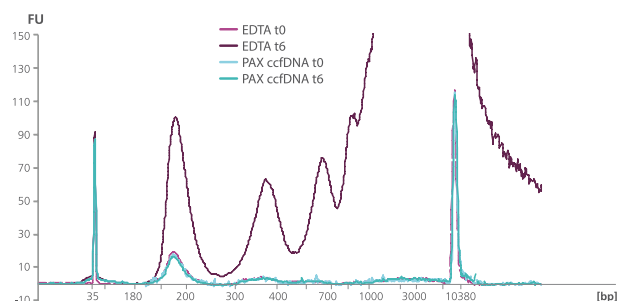

**Figure 1. PAXgene Blood ccfDNA stabilization helps prevent release of gDNA into plasma.** Whole blood was stored in EDTA tubes or PAXgene Blood ccfDNA Tubes. ccfDNA was purified from the plasma immediately following blood collection (t0) or after 6 days storage (t6) at room temperature. Eluate (1 µl) was analyzed using the Agilent High Sensitivity DNA Kit. Plasma from EDTA tubes showed an increase in apoptotic gDNA fragments, whereas plasma from PAXgene Blood ccfDNA Tubes showed a ccfDNA profile comparable to day 0.

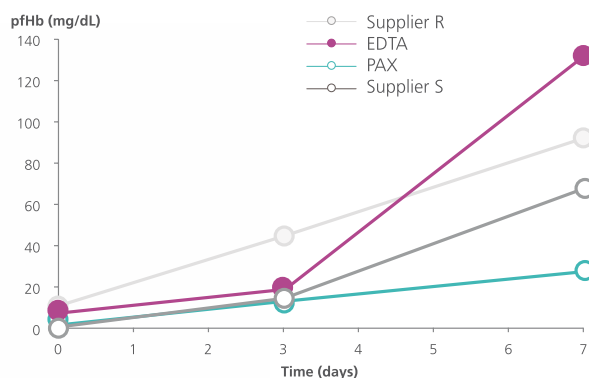

**Figure 2. PAXgene Blood ccfDNA stabilization minimizes hemolysis compared to alternative solutions.** Whole blood from 20 subjects was collected into PAXgene Blood ccfDNA Tubes, EDTA tubes and blood collection tubes designed for ccfDNA stabilization from two other suppliers. Plasma was separated directly after blood draw, after 3 days and after 7 days storage at room temperature (15–25°C). Hemolysis was assessed. Increase in sample hemolysis during sample storage at room temperature was minimized in PAXgene Blood ccfDNA Tubes compared to the other blood collection tubes.



Manual and automated solutions


# Get the most from your precious sample

## ccfDNA preparation

Choose from our ccfDNA extraction portfolio with unmatched flexibility, for high yield, high concentration ccfDNA extraction: **QIAamp MinElute ccfDNA** and **EZ1 ccfDNA kits** with scalable input volumes and workflows that can be automated on QIAcube and EZ1 Advanced XL.

**QIAasympyony SP** automates your nucleic acid preparations and seamlessly integrates into your daily work. It standardizes your workflow and delivers reproducible, high-quality data. The QIAasympyony SP enables sample preparation of nucleic acids – including DNA, ccfDNA, RNA and bacterial or viral nucleic acids – from a wide range of starting materials, such as whole blood, buffy coat, serum, plasma, urine and tissue.

### ccfDNA purification for research applications

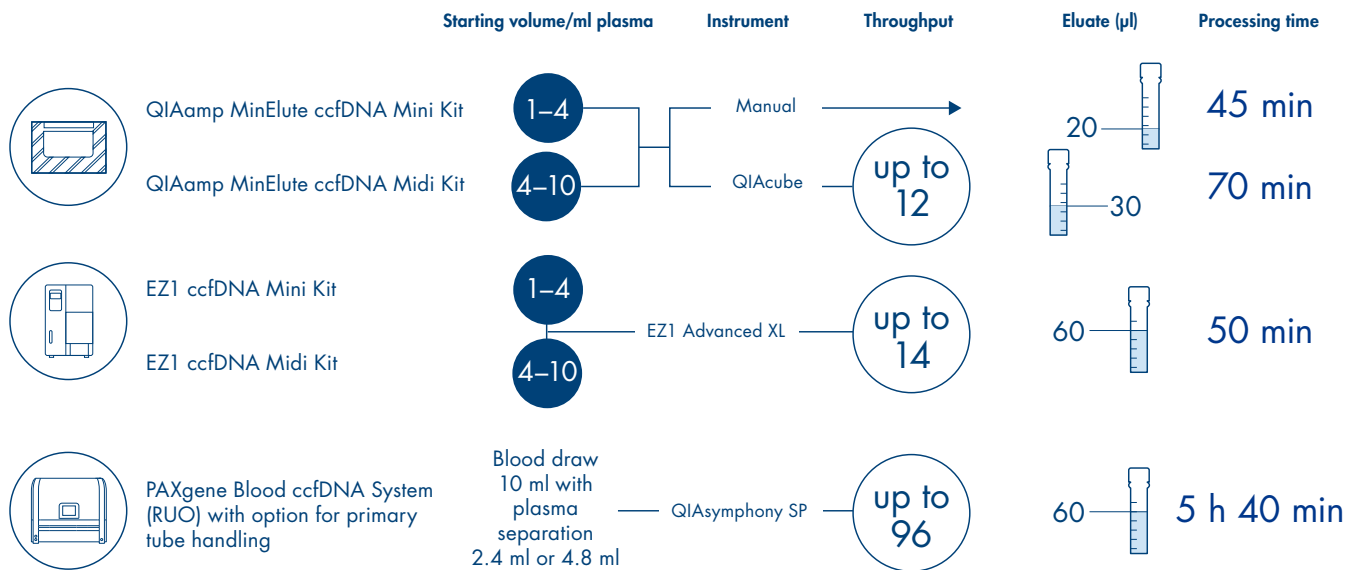

### ccfDNA purification for diagnostic applications

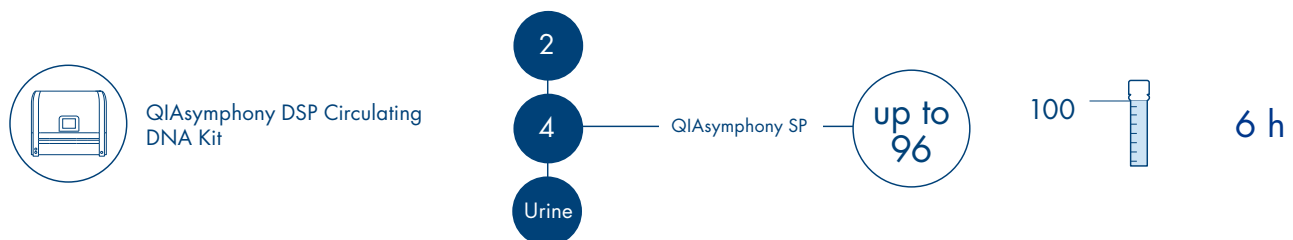


Figure 3. ccfDNA purification selection guide.

## Preparation of cell-free miRNA

The **miRNeasy Serum/Plasma Advanced Kit** efficiently purifies total RNA, including miRNA, from serum and plasma samples. The phenol-free protocol uses easy-to-automate MinElute spin column technology to allow for small elution volumes. Extracellular vesicles are lysed to allow analysis of RNA in and outside of vesicles.

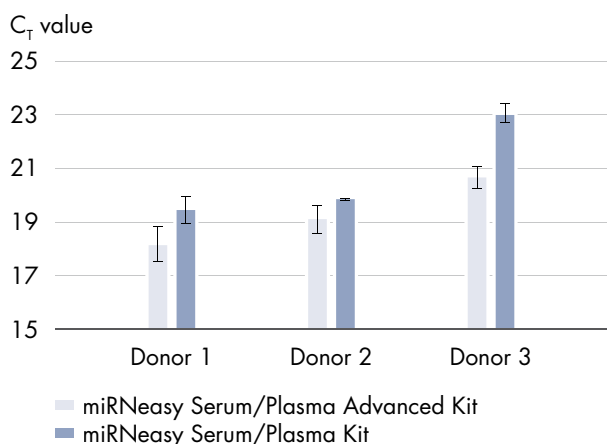


Figure 4. High recovery of miRNA without the need for phenol.

## You want them all?

The **QIAamp ccfDNA/RNA Kit** offers co-isolation of ccfDNA and ccfRNA from the same plasma sample, providing excellent yield and quality of both analytes in the same eluate. This allows you to analyze both DNA and RNA mutations in the same liquid biopsy sample or can boost detection sensitivity to discover low-frequency tumor mutations.

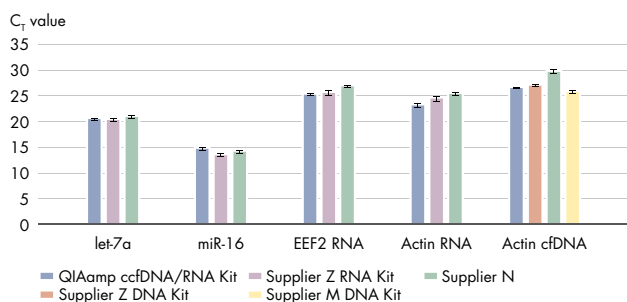


Figure 5. Competitive performance without compromising on DNA or RNA yields. Cell-free DNA/RNA was isolated in replicates from 1 ml EDTA plasma using the QIAamp ccfDNA/RNA Kit or kits from other suppliers.

## Exosome isolation kits

Explore the secrets hidden in exosomes. Enjoy spin column simplicity with ultracentrifugation quality – try the **exoRNeasy** and **exoEasy Kits**.

- Isolate extracellular vesicles in just 25 minutes
- Use high input volumes to get low abundance transcripts
- Separate vesicular vs. non-vesicular miRNA

## The CTC AdnaTest

The **AdnaTest** is a highly specific immunomagnetic cell selection system that combines an optimized antibody mixture with highly sensitive RT-PCR technology.

- Detection and molecular characterization of CTCs
- Highest CTC specificity and sensitivity

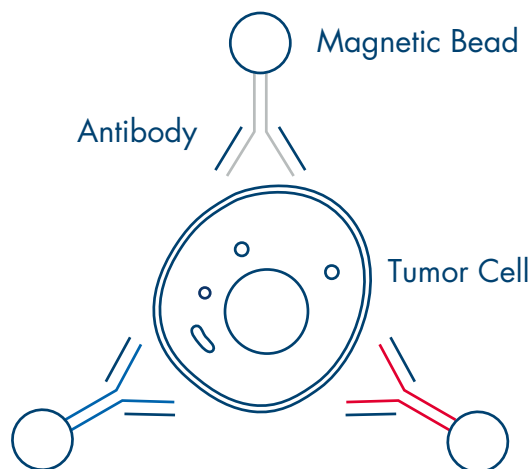


Figure 6. Tumor cell captured by 3 different antibodies ensuring highly specific isolation and detection of CTCs.

# Explore biomarkers in any liquid biopsy sample with NGS

## QIAseq Targeted Panels


Through Digital NGS powered by unique molecular indices (UMIs), **QIAseq Targeted DNA panels** allow the detection of low-frequency variants with high confidence by overcoming the issues of PCR duplicates, false positives and library bias. Target design through single primer extension (SPE) removes conventional restrictions of two-primer amplicon designs – also reducing the amount of required primers.

## QIAseq library preparation


For new biomarker discoveries or to search for aneuploidies from any cell-free DNA sample, QIAGEN's **QIAseq cfDNA Library Kits** provide highly efficient library preparation to ensure maximum sample conversion and sensitive variant detection – for whole genome or exome sequencing applications.

## QIAseq miRNA library preparation


QIAseq miRNA sequencing kits deliver high-quality data by using unique molecular indices (UMIs) and a simple gel-free workflow. **The QIAseq miRNA Library Kit** achieves greater sequencing efficiency by eliminating adapter dimers and unwanted RNA species, resulting in the most efficient use of your sequencing instrument for miRNA, piRNA and small RNA discovery.

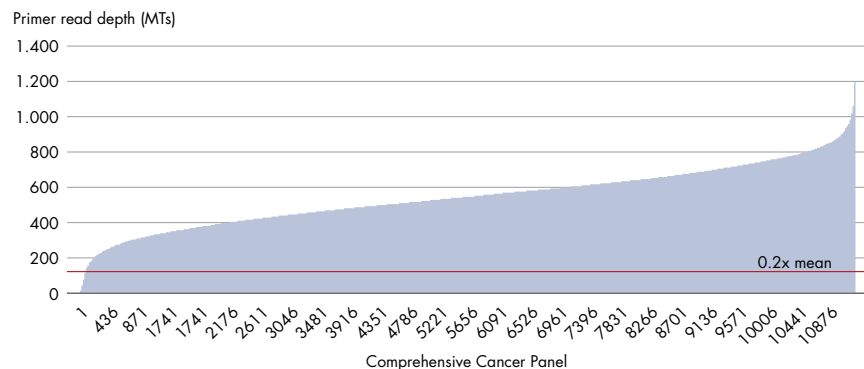

**Figure 7. Outstanding uniformity of target regions achieved by QIAseq Targeted DNA panels.** The panel achieved a uniformity of 99.5% at 0.2x of mean coverage.

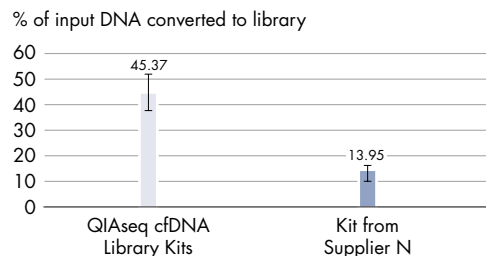


**Figure 8. Superior conversion rate of cfDNA to NGS library.** The average calculated conversion rate of the replicate samples is displayed. QIAseq cfDNA Library Kits shows significantly higher conversion rates.

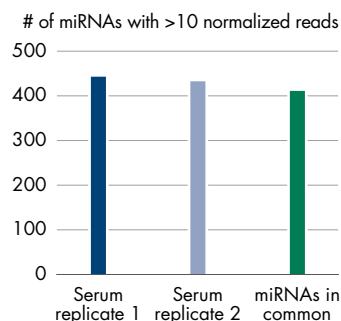

**Figure 9. The QIAseq miRNA Sequencing Kit has been designed to enhance yields from biofluids such as serum.** The graph shows robust detection of miRNA from serum samples.

PCR-based  
RNA biomarker  
profiling

# Discover disease-specific biomarker signatures

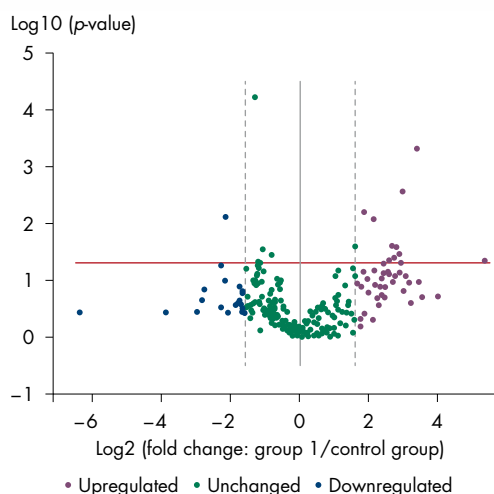

**Figure 10. Differential detection of miRNA expression using the miScript PCR System.** Volcano plot of miRNA expression changes in NSCLC patient serum samples, compared with healthy donor serum samples. Y: p-value; X: log<sub>2</sub> (fold change).

Simultaneously profile mRNA, miRNA and lncRNA

Total RNA discovery. The study of RNA has evolved from the simplicity of the central dogma of molecular biology. There are now multiple known noncoding RNA species that directly regulate gene expression. To truly understand gene expression, exploring regulatory RNA, such as miRNA and lncRNA, is key.

- miRCURY PCR Panels
- RT<sup>2</sup> PCR Arrays
- RT<sup>2</sup> lncRNA PCR Arrays


## Solutions for mRNA quantification

The QuantiNova family of qPCR kits combines top-notch performance with various in-process control features to ensure unbiased, reproducible results – even with minor changes in transcript levels.


Process and analysis errors are reduced with:

- Internal Control RNA
- Visual pipetting control
- gDNA removal
- Room temperature setup

While the **QuantiNova Multiplex PCR and RT-PCR Kits** provide the maximum information from a single PCR run, a complete kit portfolio for 1- or 2-step RT-PCR covers any qPCR need you may have for screening or confirmation – delivering instant success at the first attempt.

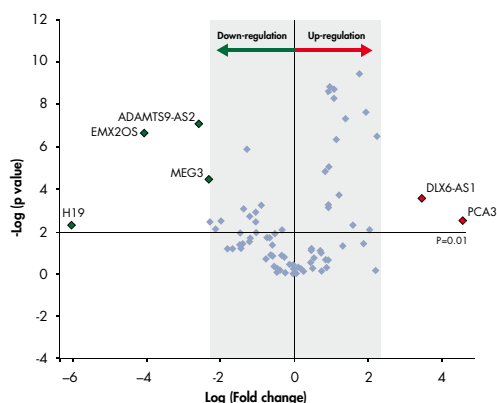

**Figure 11. Successful detection of lncRNAs in serum samples with preamplification.** Volcano plot of lncRNA gene expression changes in NSCLC serum samples compared with healthy donor serum samples. Y: p-value; X: Log<sub>2</sub> (fold change).

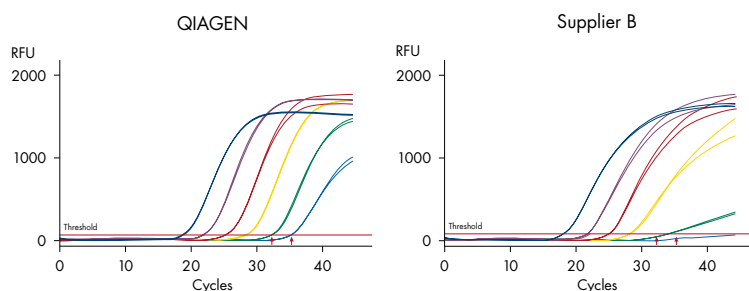

**Figure 12. Superior sensitivity of the QuantiNova Probe RT-PCR Kit.** The performance was compared to a probe RT PCR kit from supplier B using 10-fold dilutions of HeLa total RNA from 100 ng to 10 pg. Sensitivity was clearly differing for 100 and 10 pg of total RNA (red arrows).




# Concentrate on insights, not data

Data analysis  
and  
interpretation

## Variant identification and interpretation of ccfDNA

There are many advantages to using cell-free DNA for identifying new variant biomarkers and for monitoring cancer patients after drug treatment. However, data analysis can be difficult, as it requires identifying low-frequency variants at levels of 1% and lower. In addition, data interpretation for detecting potential cancer driver variants comes with its own challenges, such as identifying variants having a functional impact.

QIAGEN Bioinformatics Advanced Testing Solution, featuring **Biomedical Genomics Workbench** and **Ingenuity Variant Analysis**, combines state-of-the-art data analysis and data interpretation into a single offering. Optimized and streamlined end-to-end workflows enable high sensitivity in detection, filtering and interpretation of known and potential new cancer driver variants down to 1% (or lower) allelic fractions.

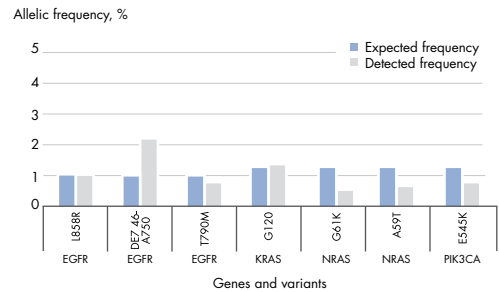

**Figure 13. Sensitive mutation detection from cfdNA.** Pathogenic variants at a ~1% allelic frequency were detected even at a moderate average read coverage of 500.

## RNA sequencing for exosomes

RNA and small RNA sequencing are powerful tools in gene and transcript expression profiling for determining underlying disease mechanisms and discovering disease-specific biomarker signatures. The RNA-Seq Explorer Solution enables you to start from FASTQ and align, analyze, QC, normalize and determine differential expression. Now you can identify the biology associated with differentially expressed genes and isoforms of varying expressions without being a bioinformatician.

## Understand complex 'omics data with IPA

**Ingenuity Pathway Analysis (IPA)** is a powerful analysis and search tool that uncovers the significance of 'omics data and identifies new targets or candidate biomarkers within the context of biological systems. IPA has broadly been adopted by the life science research community and is cited in thousands of articles for the analysis, integration and interpretation of data derived from 'omics experiments, such as RNA-seq, small RNA-seq, microarrays including miRNA and SNP, metabolomics, proteomics and small-scale experiments.



## Ordering Information

Product	Contents	Cat. no.
<b>Blood collection and stabilization</b>		
PAXgene Blood ccfDNA Tubes (100)	100 blood collection tubes (10 ml)	768115
<b>Preparation</b>		
QIAamp MinElute ccfDNA Mini or Midi Kit (50)	Purification of free-circulating DNA from 1–4 ml (Mini) or 4–10 ml (Midi) serum or plasma, automatable on QIAcube	55204 55284
EZ1 ccfDNA Mini or Midi Kit (48)	Automated purification of circulating cell-free DNA from 1–4 ml (Mini) or 4–10 ml (Midi) serum or plasma on the EZ1 Advanced XL System	954134 954154
QIAasymphony PAXgene Blood ccfDNA Kit (192)	Purification of ccfDNA from plasma using the QIAasymphony SP with optimized chemistry and input volumes for PAXgene Blood ccfDNA Tubes	768536
miRNeasy Serum/Plasma Advanced Kit (50)	Phenol-free isolation of total RNA including miRNA from serum and plasma samples	217204
QIAamp ccfDNA/RNA Kit (50)	Co-purification of cell-free DNA and RNA from plasma and serum	55184
exoEasy Maxi Kit (20)	Purification of exosomes and other extracellular vesicles (EVs) from plasma, serum and cell culture supernatant	76064
exoRNeasy Serum/Plasma Midi or Maxi Kit (50)	Purification of RNA from exosomes and other extracellular vesicles from up to 1 ml (Midi) or 4 ml (Maxi) serum or plasma	77044 77064
miRCURY Exosome Serum/Plasma Kit	Enrichment of exosomes and other extracellular vesicles from serum or plasma	76603
miRCURY Exosome Cell/Urine/CSF Kit	Enrichment of exosomes and other extracellular vesicles from biofluids such as cerebrospinal fluid, urine or cell culture media	76743
EZ1 Advanced XL, System	Robotic workstation for automated purification of nucleic acids from up to 14 samples using EZ1 Kits	9001874
QIAcube	Robotic workstation for automated purification of DNA, RNA, or proteins using QIAGEN spin-column kits	9001293
QIAasymphony SP	For fully integrated automated sample preparation of 1–96 samples	9001297
<b>NGS analysis</b>		
QIAseq Targeted DNA Panels	For confident detection of low-frequency variants using unique molecular indices	Varies
QIAseq cfDNA Library T Kit	NGS library preparation from ccfDNA for whole genome or exome sequencing on Ion Torrent sequencers	1102308
QIAseq cfDNA Library Kit	NGS library preparation from ccfDNA for whole genome or exome sequencing on Illumina® sequencers	180015

## Ordering Information

Product	Contents	Cat. no.
<b>NGS analysis</b> QIAseq miRNA Library Kits*	For gel-free miRNA sequencing library preparation	331502
<b>PCR-based analysis</b> miRCURY PCR Panels	SYBR Green-based, real-time PCR profiling of miRNAs using the miScript® PCR System	Varies
RT <sup>2</sup> Profiler PCR Arrays	Pathway-focused gene expression analysis using laboratory-verified assays	Varies
RT <sup>2</sup> lncRNA Profiler Arrays	lncRNA expression analysis using laboratory-verified SYBR® Green qPCR assays	Varies
QuantiNova PCR Kits	Highly sensitive, specific, and ultrafast real-time PCR	Varies
<b>Data analysis</b> Biomedical Genomics Workbench	For analyzing, comparing and visualizing human hereditary and cancer NGS data for biomarker discovery	832070
CLC Genomics Workbench	For analyzing, comparing, and visualizing NGS data	832000
<b>Biological interpretation</b> Ingenuity Variant Analysis	For identifying disease-causing variants from human whole genome, exome and gene panel next-generation sequencing studies	Varies
Ingenuity Pathway Analysis (IPA)	For modeling, analyzing, and understanding complex 'omics data	Varies

\*Various options are available; please inquire.

Product	Contents	Cat. no.
<b>ccfDNA preparation for diagnostic applications</b> QIAasymphony DSP Circulating DNA Kit	For fully automated purification of human circulating cell-free DNA from human plasma and urine using the QIAasymphony SP instrument	937556

### Visit [qiagen.com](http://qiagen.com) today to find out more about Liquid Biopsy products!

The QIAasymphony DSP Circulating DNA Kit is for in vitro diagnostic use. All other products listed in this brochure are intended for molecular biology applications. These products are not intended for the diagnosis, prevention or treatment of a disease.

For up-to-date licensing information and product-specific disclaimers, see the respective QIAGEN kit handbook or user manual. QIAGEN kit handbooks and user manuals are available at [www.qiagen.com](http://www.qiagen.com) or can be requested from QIAGEN Technical Services or your local distributor.

---

Trademarks: QIAGEN®, Sample to Insight®, QIAamp®, QIAcube®, QIAseq™, QIASymphony®, EZ1®, Ingenuity®, MinElute®, miScript®, QuantiFast®, QuantiNova®, Rotor-Gene®, RT<sup>2</sup>™, RT<sup>2</sup> Profiler™ (QIAGEN Group); PAXgene® (PreAnalytiX); Illumina® (Illumina, Inc.); SYBR® (Thermo Fisher Scientific Inc.). Registered names, trademarks, etc. used in this document, even when not specifically marked as such, are not to be considered unprotected by law.

1107666 PROM-10170-002 09/2017 © 2017 QIAGEN, all rights reserved.

Ordering [www.qiagen.com/shop](http://www.qiagen.com/shop) | Technical Support [support.qiagen.com](http://support.qiagen.com) | Website [www.qiagen.com](http://www.qiagen.com)